



Human recombinant Endocan (50 kDa)

Carrying a unique DS / DS chain

Essential Notes

Cat. Number : LIP-1001

Biomolecule : endocan

Concentration: 100 µg/mL

Size :

10 µg – LIP-1001-10U

50 µg – LIP-1001-50U

Formulation : PBS pH 7.4

Storage : -20°C

Applications : ELISA, Bioassays, standard for quantification

FOR RESEARCH USE ONLY

Description

Endocan also called endothelial cell specific molecule-1 (ESM-1) is a proteoglycan secreted by endothelial cells, and up-regulated by proangiogenic factors (VEGF, FGF-2) and pro-inflammatory cytokines (TNF- α , IL-1). Endocan plays a role in the angiogenesis and neovessel formation in association with VEGF, in the control of inflammation by inhibiting the leukocyte diapedesis, and in the regulation of the cellular proliferation by regulating growth factor activities through its glycan.

Source

The recombinant human endocan / ESM-1 is produced by a HEK293 cell line overexpressing human endocan, obtained by transformation of HEK293 cells with cDNA encoding for human endocan. Recombinant endocan was purified by ion exchange chromatography followed by affinity chromatography using a monoclonal anti-human endocan antibody.

Molecular Mass

As a result of post-translational modification, human endocan / ESM-1 is secreted as a 50 kDa proteoglycan.

Formulation

The purified recombinant endocan / ESM-1 is diluted in a solution of phosphate buffer saline (PBS) at 100 µg/mL.

Storage


Upon reception, store at -20 to -80°C for long term storage. Avoid repeated freeze-thaw cycles.

Applications

Bioassays : Optimal dilutions should be determined by each laboratory for each application.

Western Blot (WB) : Can be used as a standard for molecular weight (50 kDa).

ELISA : Can be used for the quantification of human endocan in serum, plasma or culture supernatant as a standard.

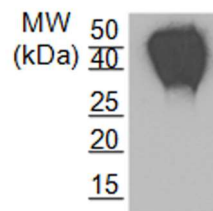

Western blot analysis of Recombinant Human Endocan. The 50 kDa endocan proteoglycan was detected using the purified mAb MEP14 (LIA-1001).

Address : BIOTHELIS SAS

Pavé du moulin, 59260 Hellemmes-Lille, France

Supports & Orders: contact@biothelis.fr - Tel: (33) 374 098 262

Last Updated Version: Rev 00 – January 2021

■ Bibliography related to endocan

- Gaudet A, et al. Cleaved endocan acts as a biologic competitor of endocan in the control of ICAM-1-dependent leukocyte diapedesis. *J Leukoc Biol.* 2020 May;107(5):833-841.
- Gaudet A, et al. Decrease of the plasmatic endocan cleavage ratio is associated with the hyperinflammatory phenotype of acute respiratory distress syndrome. *Crit Care.* 2019 Jul 11;23(1):252.
- Gaudet A, et al. Endocan regulates acute lung inflammation through control of leukocyte diapedesis. *J Appl Physiol* (1985). 2019 Sep 1;127(3):668-678.
- Ying J, Zhou D, Gu T, Huang J. Endocan, a Risk Factor for Developing Acute Respiratory Distress Syndrome among Severe Pneumonia Patients. *Can Respir J.* 2019 Apr 1;2019:2476845.
- Gaudet A, et al. Parmentier E, De Freitas Caires N, Portier L, Dubucquoi S, Poissy J, Duburcq T, Hureau M, Lassalle P, Mathieu D. Impact of acute renal failure on plasmatic levels of cleaved endocan. *Crit Care.* 2019 Feb 19;23(1):55
- De Freitas Caires N, et al. Endocan, sepsis, pneumonia, and acute respiratory distress syndrome. *Crit Care.* 2018 Oct 26;22(1):280.
- Yassine H, et al. The non-glycanated endocan polypeptide slows tumor growth by inducing stromal inflammatory reaction. *Oncotarget.* 2015 Feb 20;6(5):2725-35.
- De Freitas Caires N, et al. Identification of a 14 kDa endocan fragment generated by cathepsin G, a novel circulating biomarker in patients with sepsis. *J Pharm Biomed Anal.* 2013 May 5;78-79:45-51.
- Depontieu F, et al. Development of monoclonal antibodies and ELISA specific for the mouse vascular endocan. *J Immunol Methods.* 2012 Apr 30;378(1-2):88-94.
- Sarrazin et al. (2010) Characterization and binding activity of the chondroitin/dermatan sulfate chain from Endocan, a soluble endothelial proteoglycan. *Glycobiology.* 20:1380-1388.
- Depontieu F, et al. Loss of Endocan tumorigenic properties after alternative splicing of exon 2. *BMC Cancer.* 2008 Jan 18;8:14.
- Bechard et al. (2001) Endocan is a novel CS/DS proteoglycan that promotes HGF/SF mitogenic activity. *J. Biol. Chem.* 276:48341-48349.
- Bechard D, et al. Characterization of the secreted form of endothelial-cell-specific molecule 1 by specific monoclonal antibodies. *J Vasc Res.* 2000 Sep-Oct;37(5):417-25.

■ Background

Endocan, also known as endothelial cell-specific molecule (ESM-1), was originally discovered by Lassalle and collaborators in endothelial cells. Structurally, endocan is a dermatan sulfate proteoglycan of 50 kDa that is freely circulating in blood. Endocan binds CD11a/CD18 integrin (also called LFA-1 for Leukocyte Function-associated Antigen-1) on human leukocytes inhibiting consequently its binding to ICAM-1 and transendothelial migration. Moreover, endocan has been recently described as a biomarker of tip cells and neoangiogenesis. The expression of endocan is upregulated by pro-inflammatory molecules such as tumor necrosis factor alpha, and pro-angiogenic molecules such as vascular endothelial growth factor and fibroblast growth factor 2. Endocan binds via its dermatan sulfate chain to hepatocyte growth factor/ scatter factor. Endocan appears as a pertinent biomarker of endothelial dysfunction.

■ Companion products

- Anti-human endocan/ESM-1 mAb (N-ter) ; clone MEP08 : [LIA-0901](#)
- Anti-murine endocan/ESM-1 mAb (C-ter) ; clone MEP14 : [LIA-1001](#)
- Anti-human endocan/ESM-1 mAb (N-ter) ; clone MEP21 : [LIA-0902](#)
- Human recombinant endocan/ESM-1 (20 kDa) : [LIP-1002](#)
- DIYEK H1 (Do It Yourself Elisa Kit for Human Endocan quantification) : [LIK-1101](#)
- DIYEK C1 (Do It Yourself Elisa Kit for human Endocan and Cathepsin G-cleaved Endocan) : [LIK-1501](#)
- JDIEK H1 (Just Do It Elisa Kit for Human Endocan quantification) : [LIK-1201](#)

Not intended for use as a therapeutic agents or in diagnostic procedures. Not for use in humans or animals.